



Staffan Preutz has led innovation in optics from Polaris rimless eyewear in the 1970s, to the multi-award winning SP and now the ingenious Staffan Preutz Design.

Staffan Preutz has achieved global success through an enviable reputation for quality, knowledge and vision.

"Lightweight tailored eyewear as individual as you are"

DESIGN

Personalisation

Personalise Staffan Preutz Design luxury eyewear to match your customers' style:

- STEP 1. Choose frame shape.
- STEP 2. Tailor size to fit face shape.
- STEP 3. Select frame colour.
- STEP 4. Select side design & colour.
- STEP 5. (optional) Personalise with features.

Tailored Fit

Staffan Preutz Design eyewear can be tailored to fit the wearer perfectly using measurements taken by the optician.

Material Innovation

The NXT Plastic front sheet is:

- Strong and durable but flexible
- Lightweight
- Solvent resistant
- Stable over time maintains its curvature



Pure Optics LTD, 4 Fitzhamon Court, Wolverton Mill, Milton Keynes, MK12 6BB, England Tel: +44 (0)1908 222002 Fax: +44 (0)1908 222771 office@pureoptics.co.uk www.polariseyewear.com



SPD 13 SIDE COLOURS





A minimal timeless design formed from lightweight and flexible titanium. **SPD 12** SIDE COLOURS (TITANIUM)








The SPD 11 boasts a clean style, bold shapes and beautifully curved titanium sides.


SPD 11 SIDE COLOURS (TITANIUM)




www.polariseyewear.com

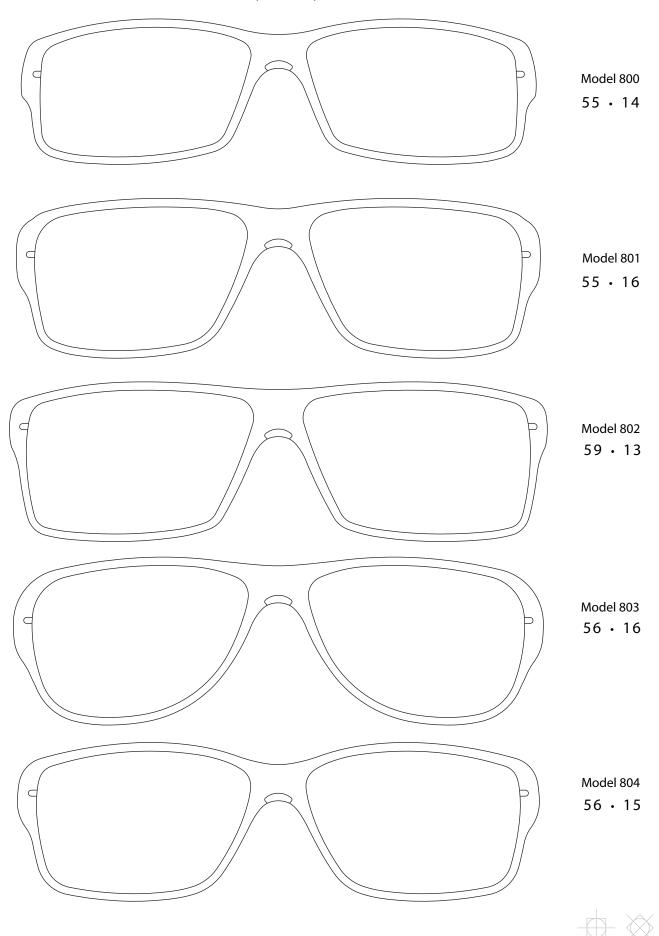




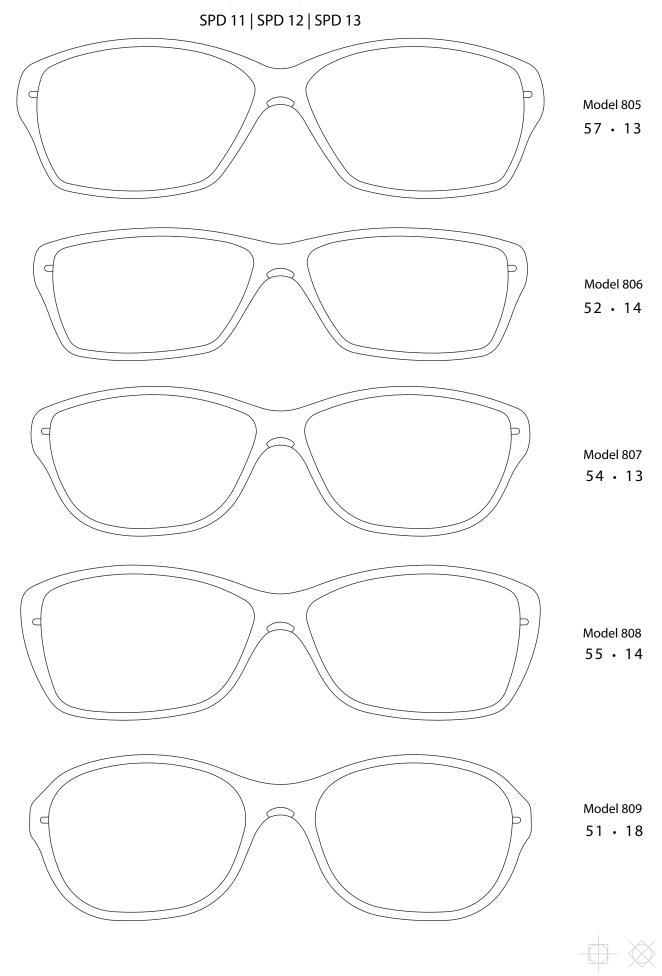

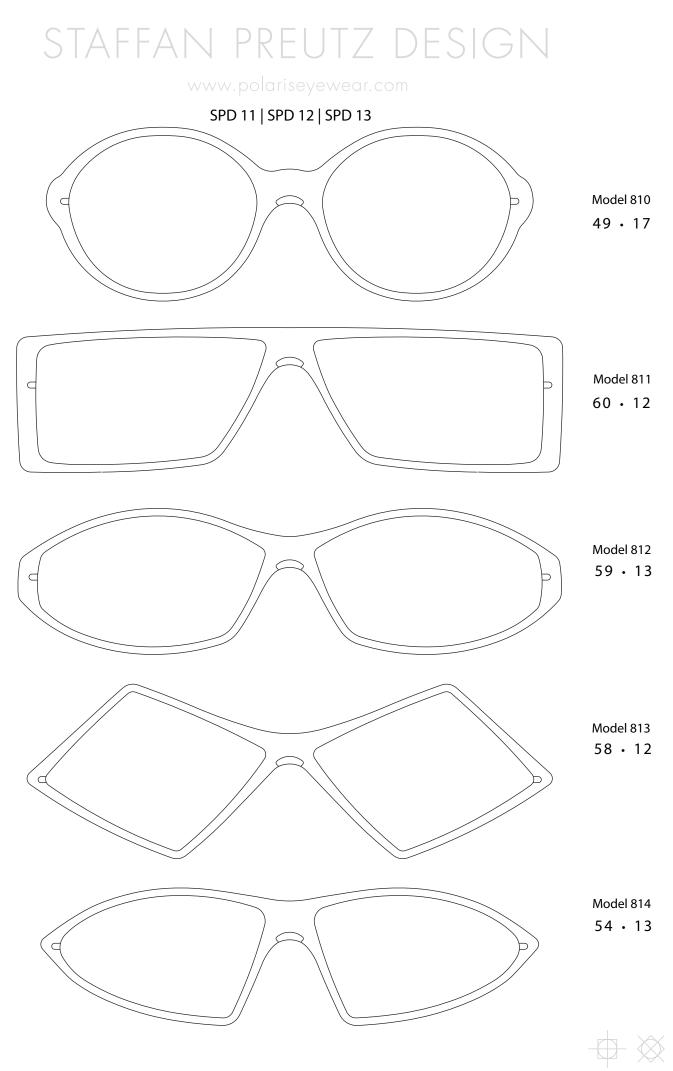



www.polariseyewear.com

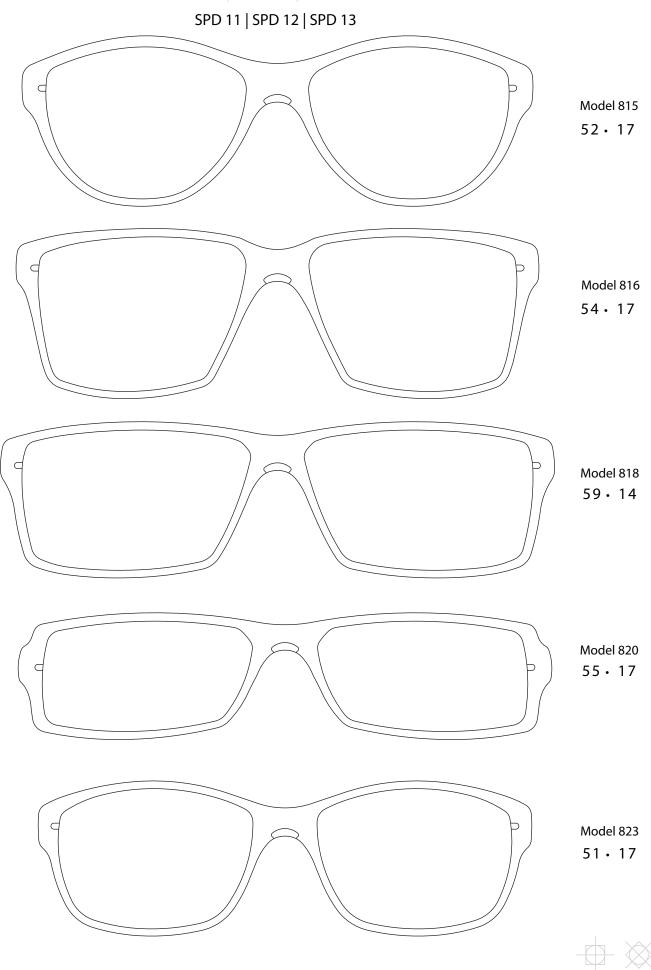


SPD11V0216

www.polariseyewear.com

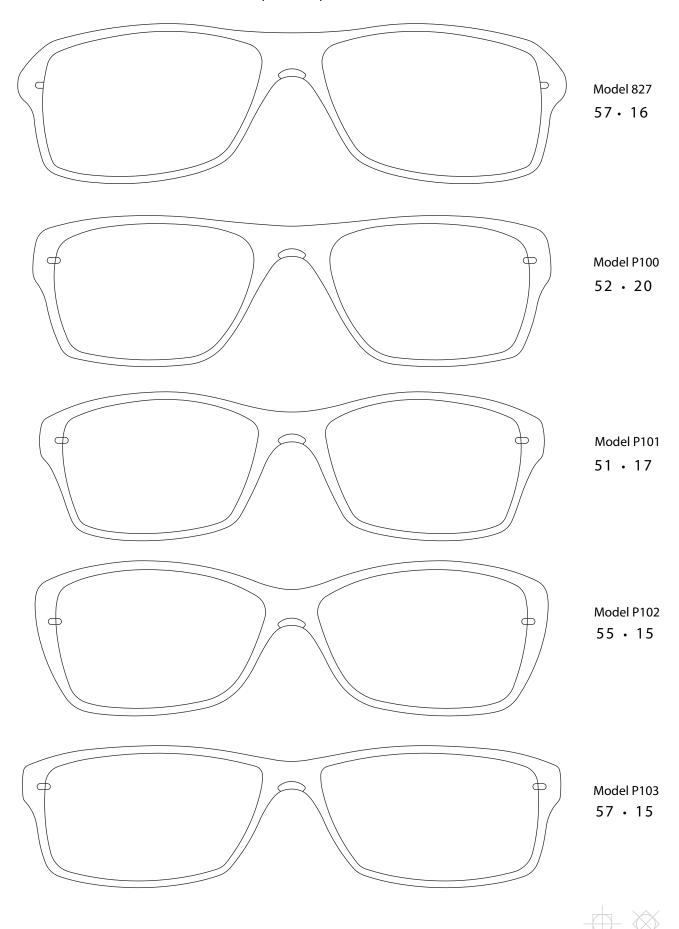


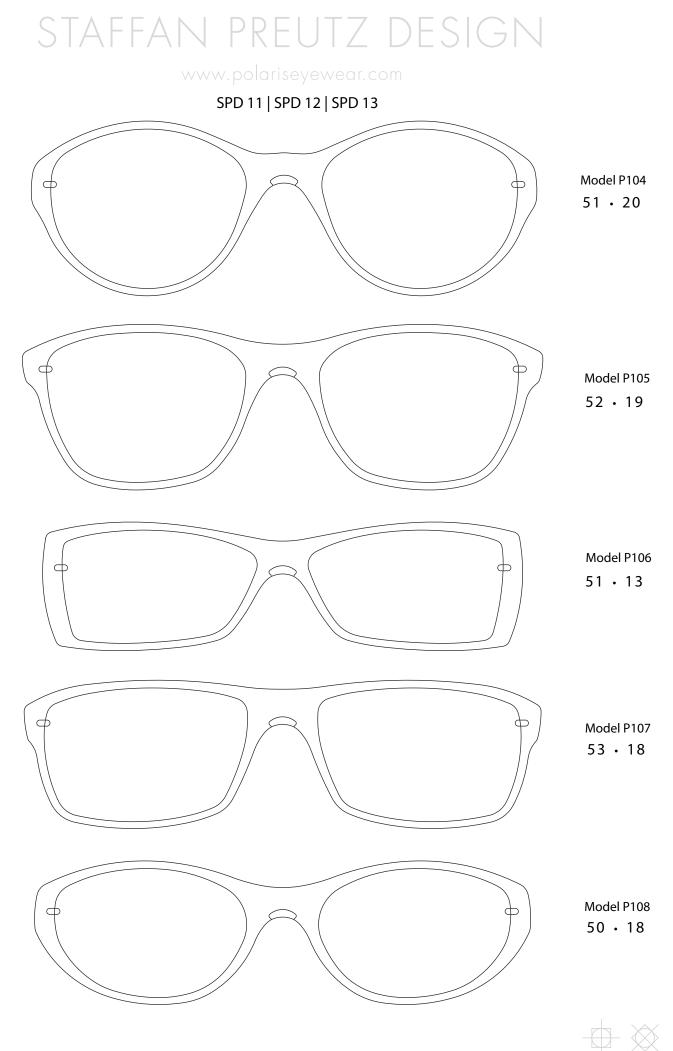

www.polariseyewear.com

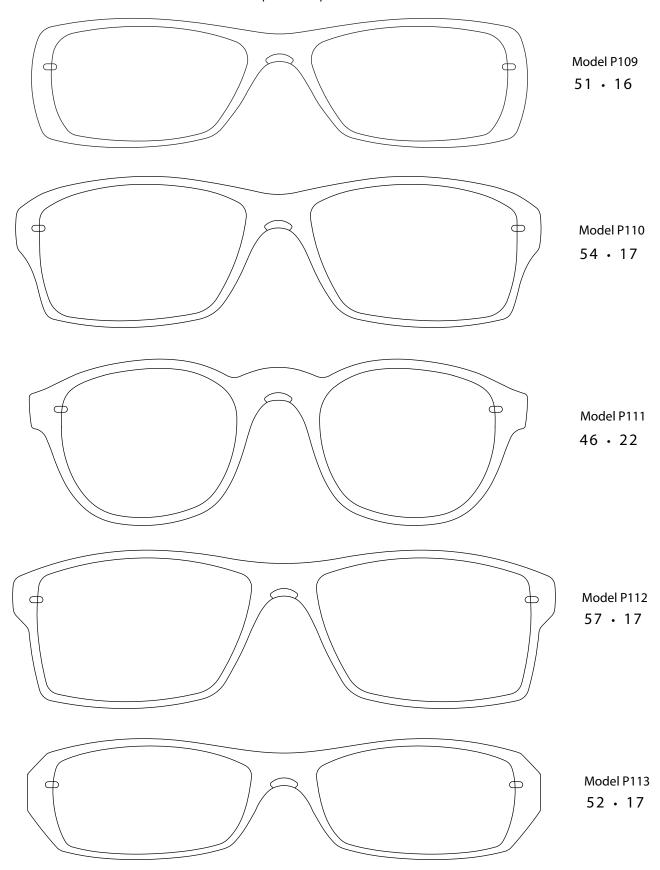


www.polariseyewear.com

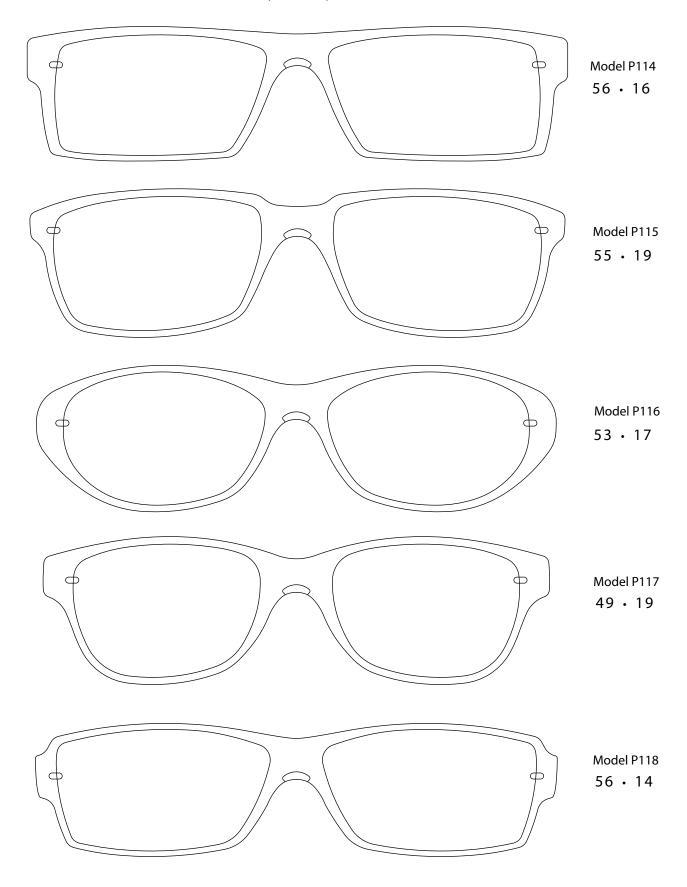





www.polariseyewear.com

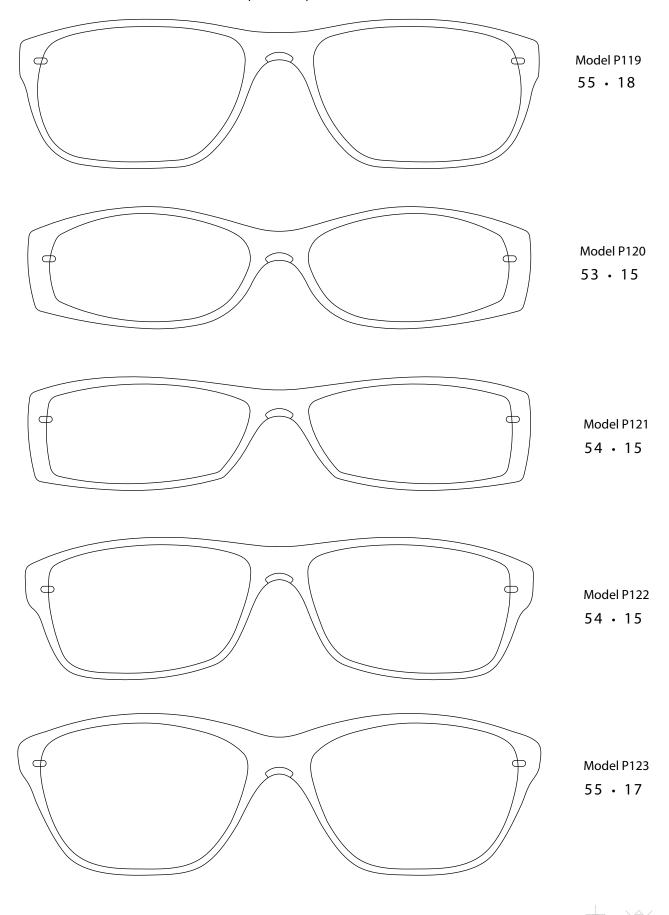


www.polariseyewear.com

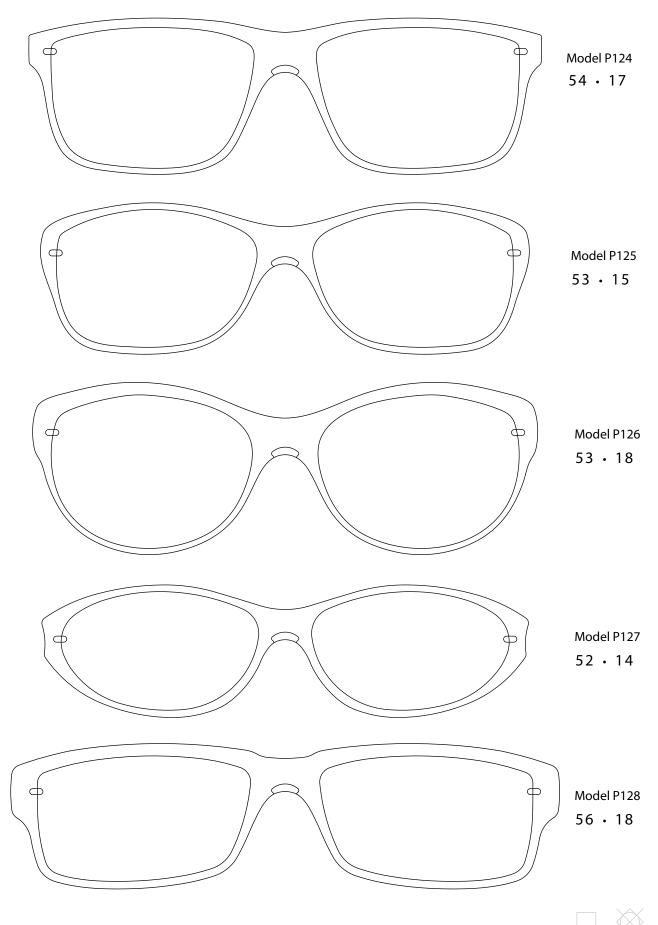


www.polariseyewear.com

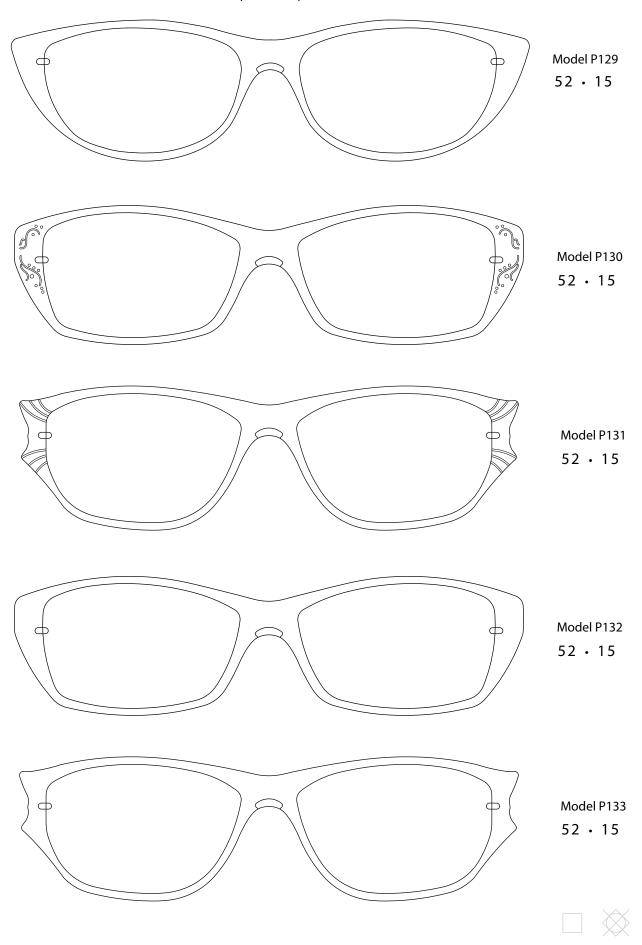


www.polariseyewear.com

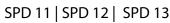

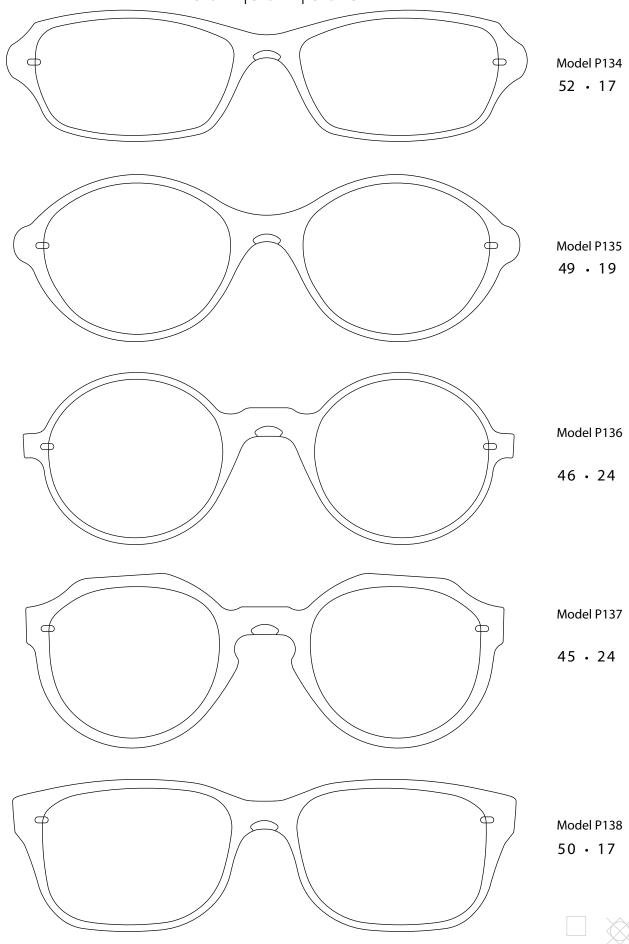

SPD13V0616

www.polariseyewear.com

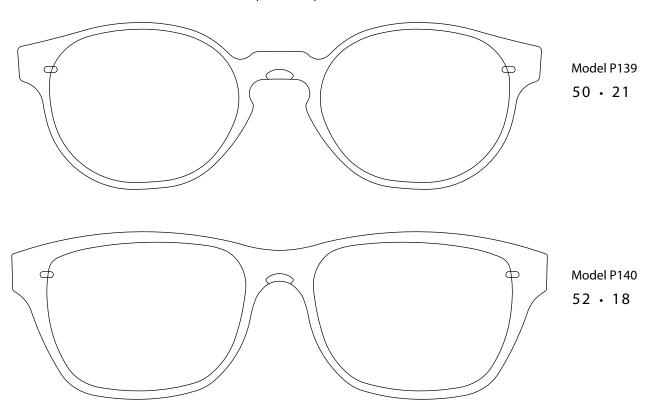





A feminine collection that adds elegance to the light and stylish design.



SPD 10 SIDE COLOURS (TITANIUM)



A feminine collection with a strong design and contrasting polished and matt finishes.

SPD 8 SIDE COLOURS (TITANIUM)









A contemporary design with a variety of sides to choose from.







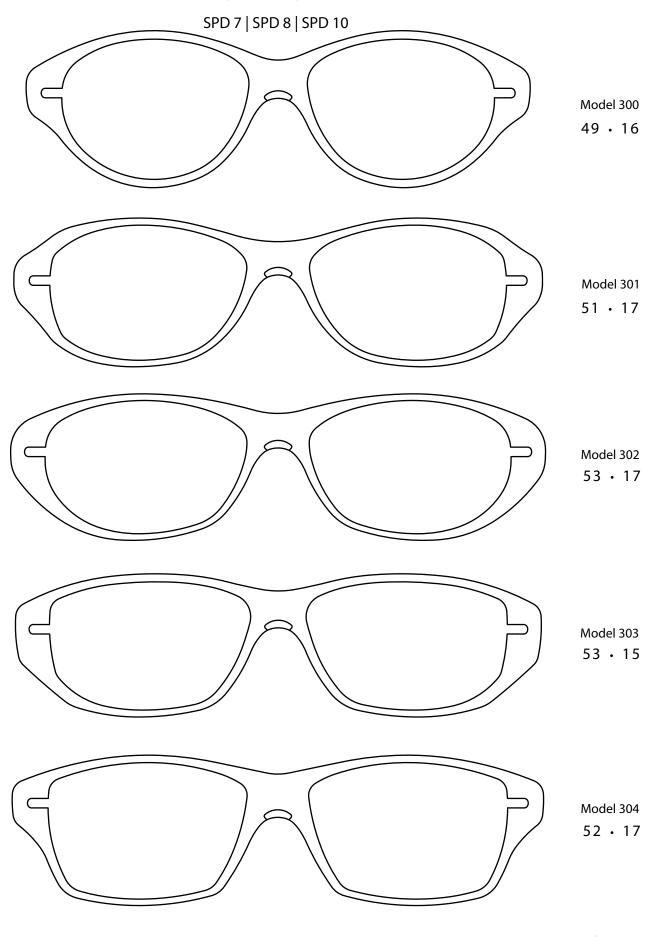



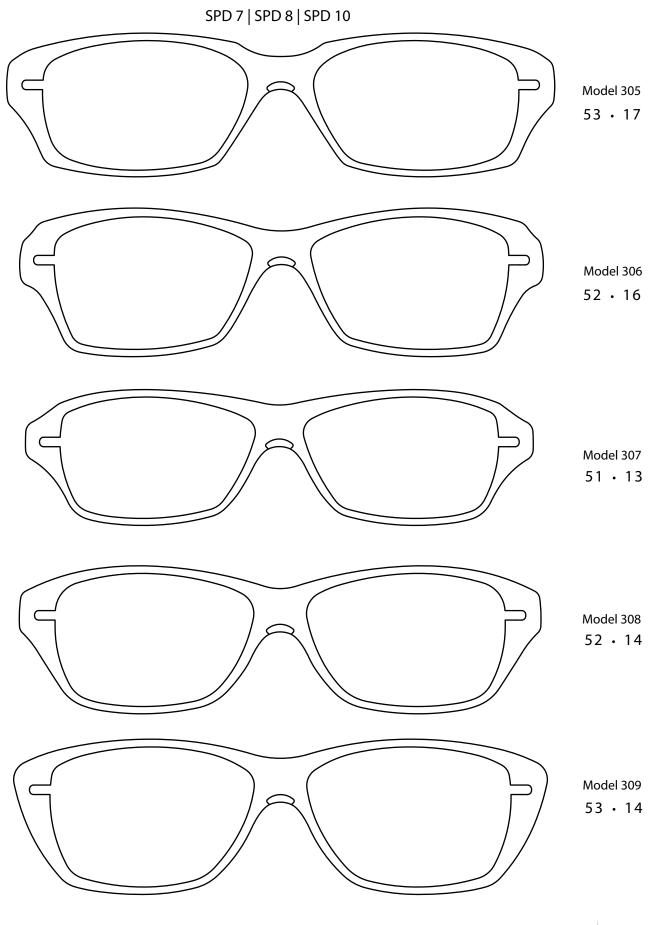

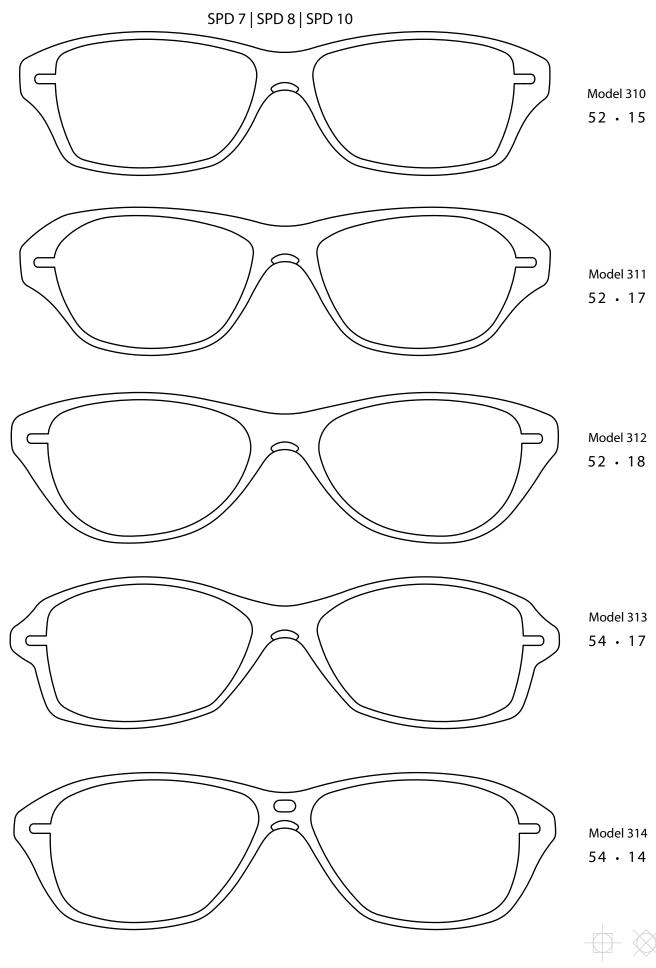


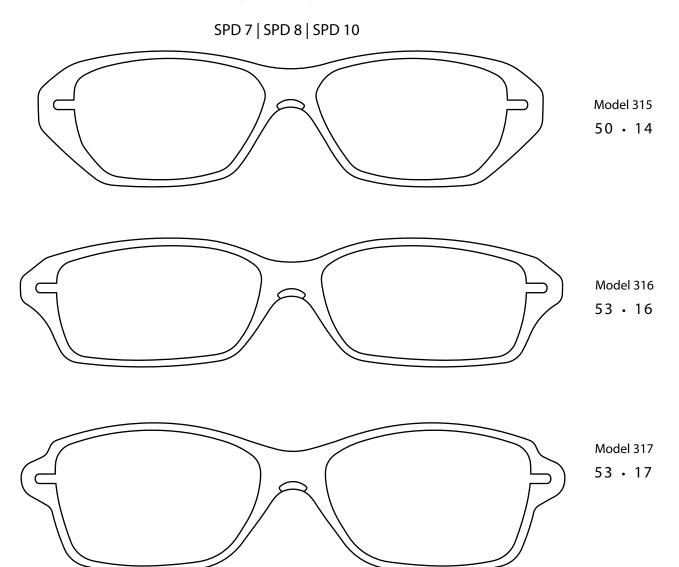




www.polariseyewear.com


SPD 7 | SPD 8 | SPD 10

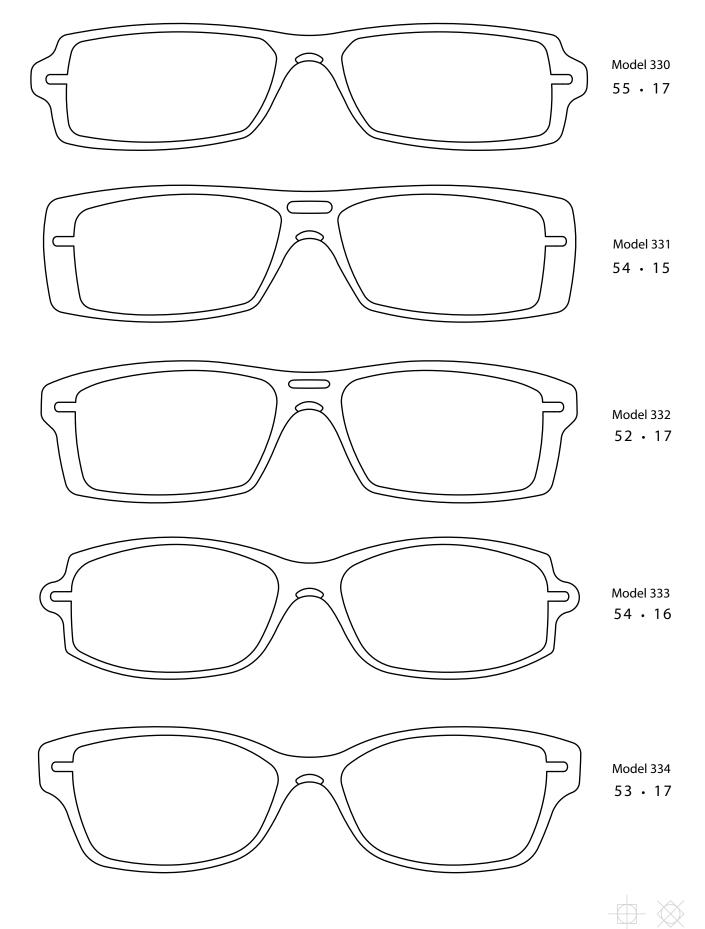


SPD7V0216

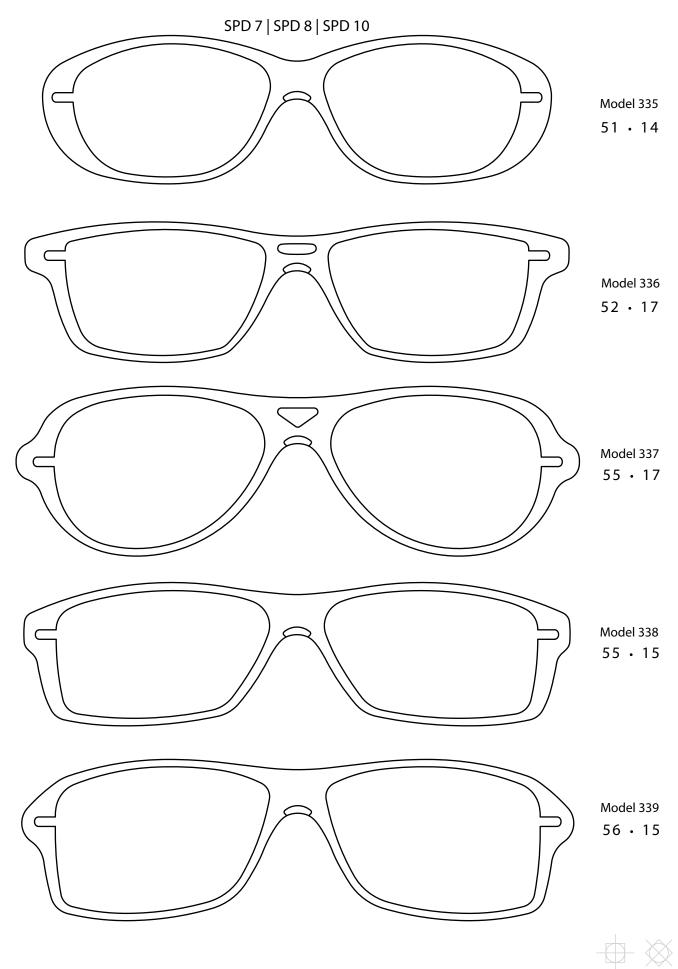


www.polariseyewear.com

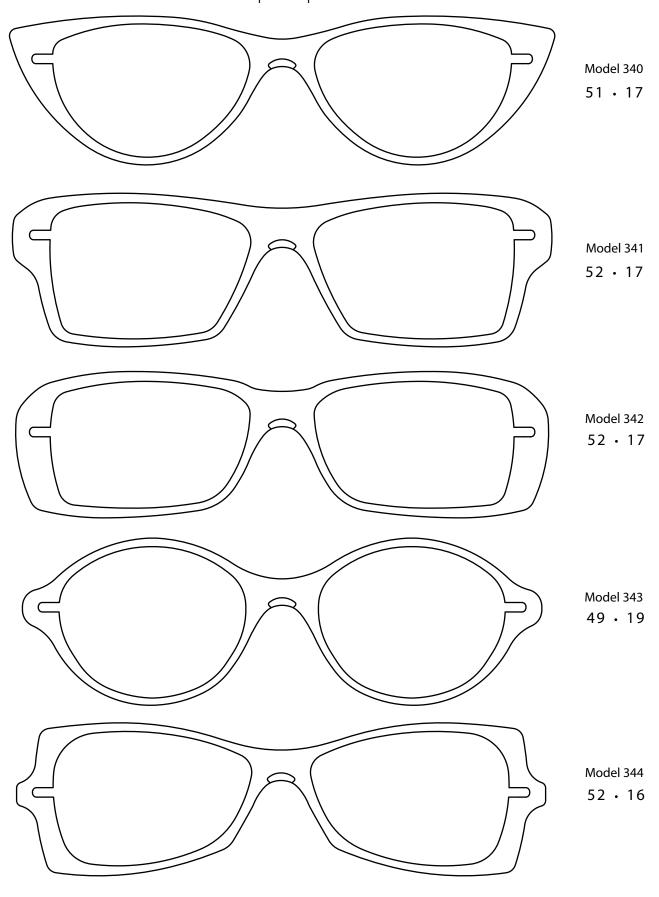


www.polariseyewear.com

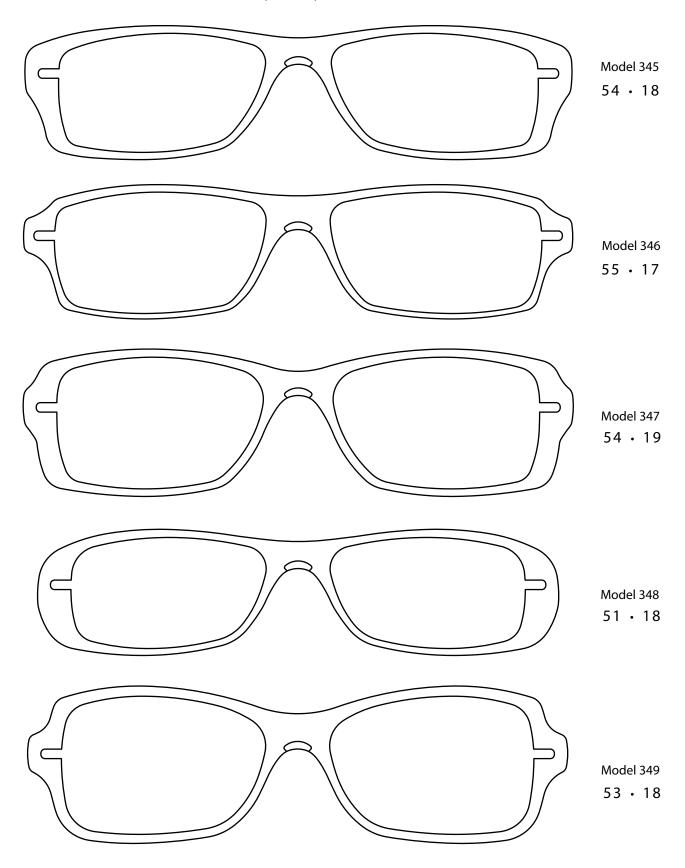



www.polariseyewear.com

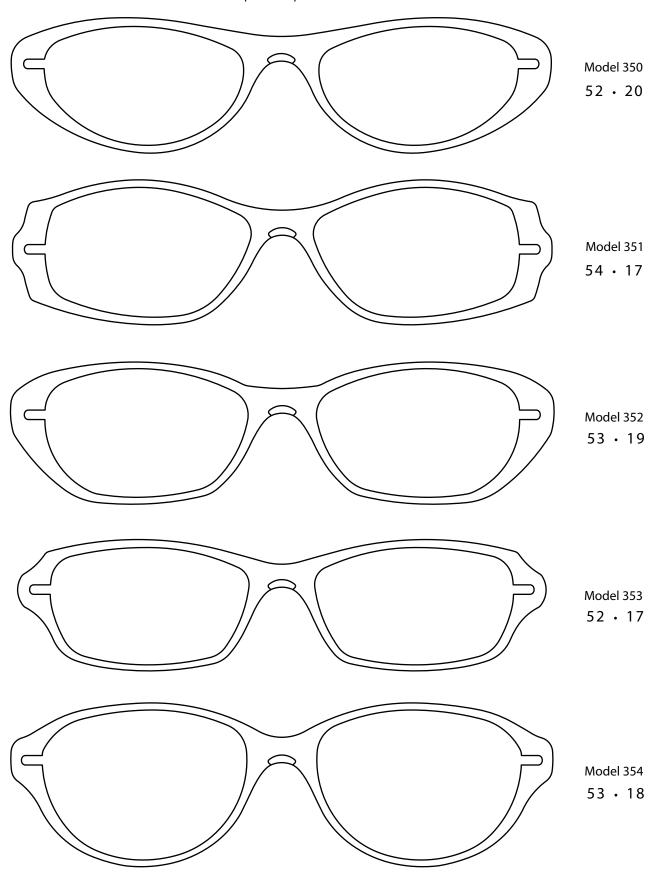

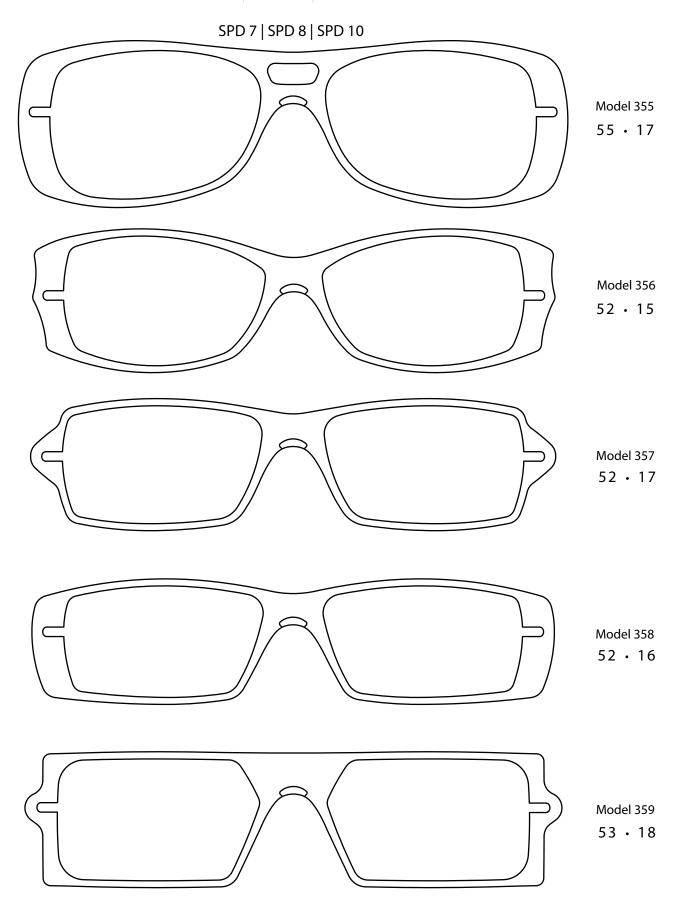


www.polariseyewear.com

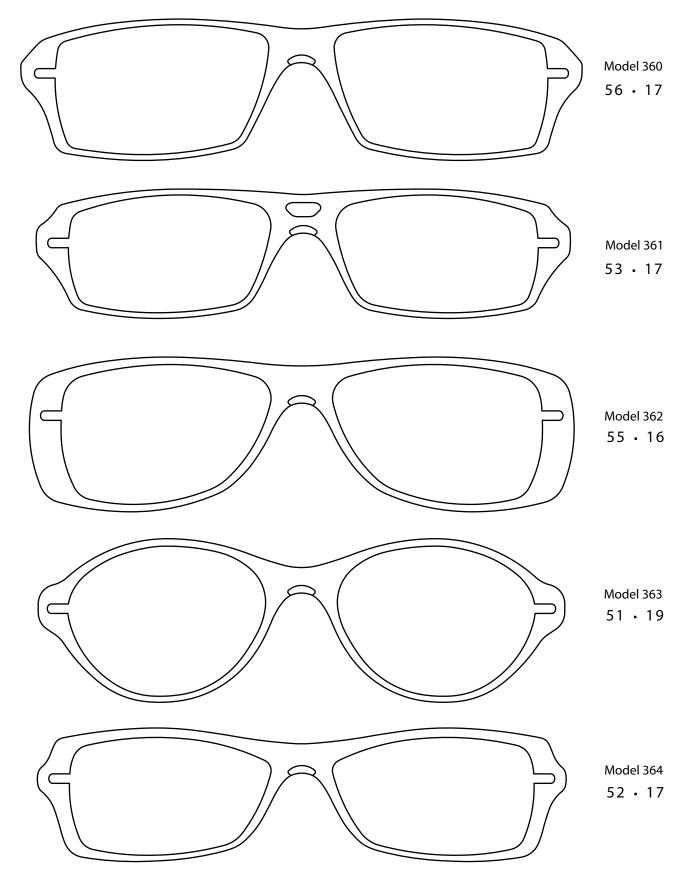

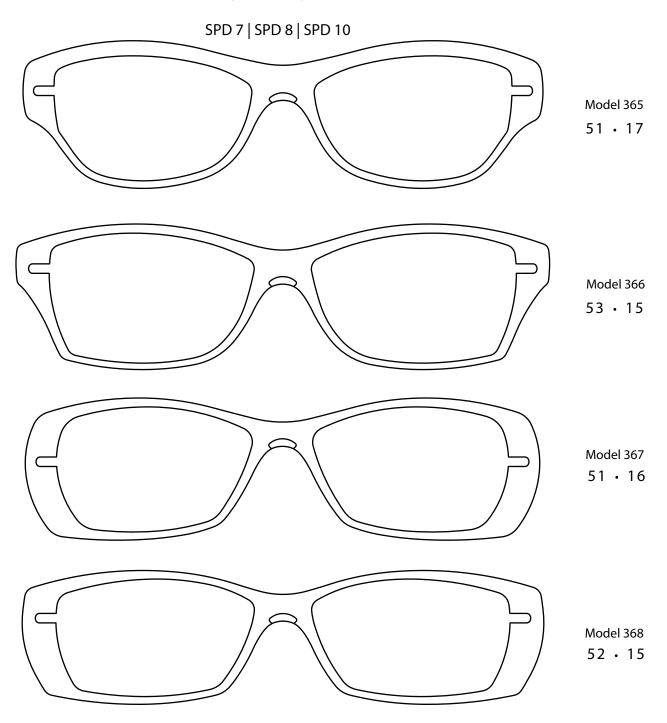





www.polariseyewear.com

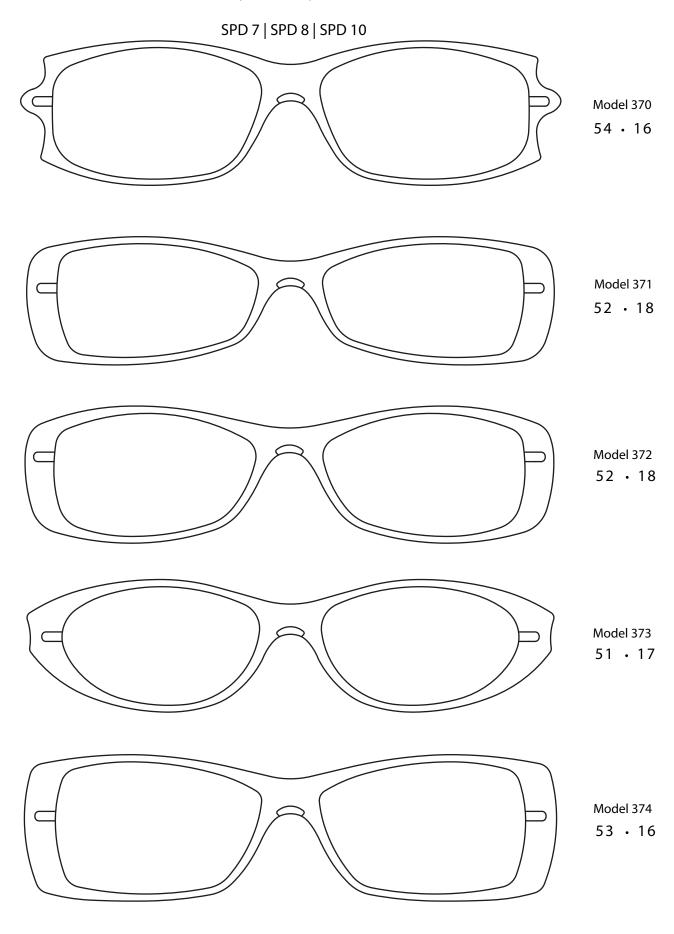

370V5014

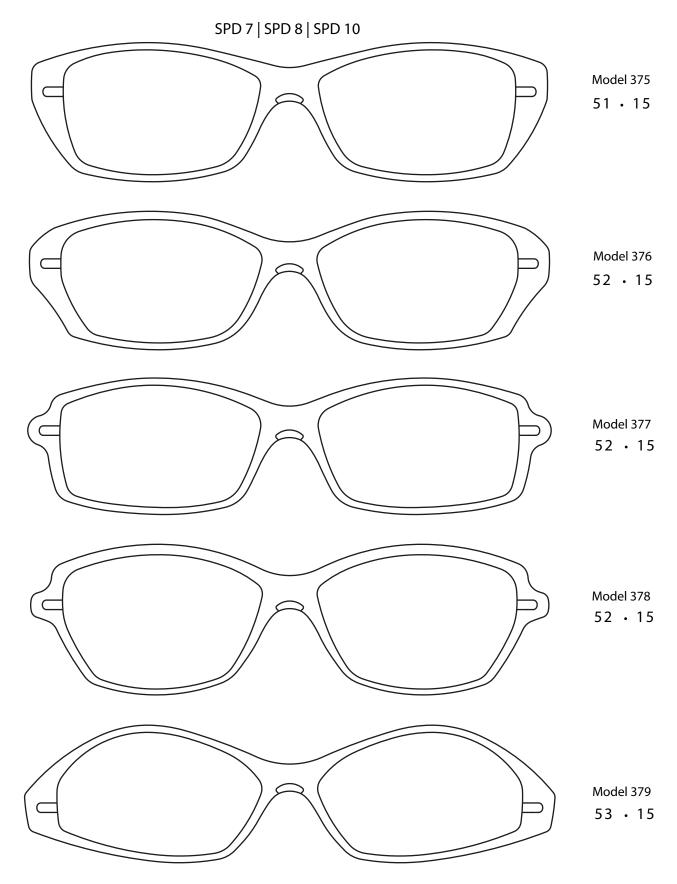



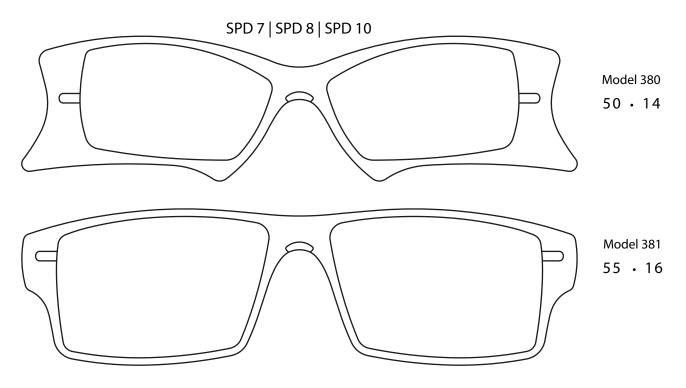



SPD 9




A collection that allows ultimate personalisation with the ability to choose both transparent and bold colours as well as different side designs.




NOTE: Due to the printing process, colours reproduced in this catalogue may vary from the actual product. Products are not shown to actual size.


www.polariseyewear.com

SPD 9

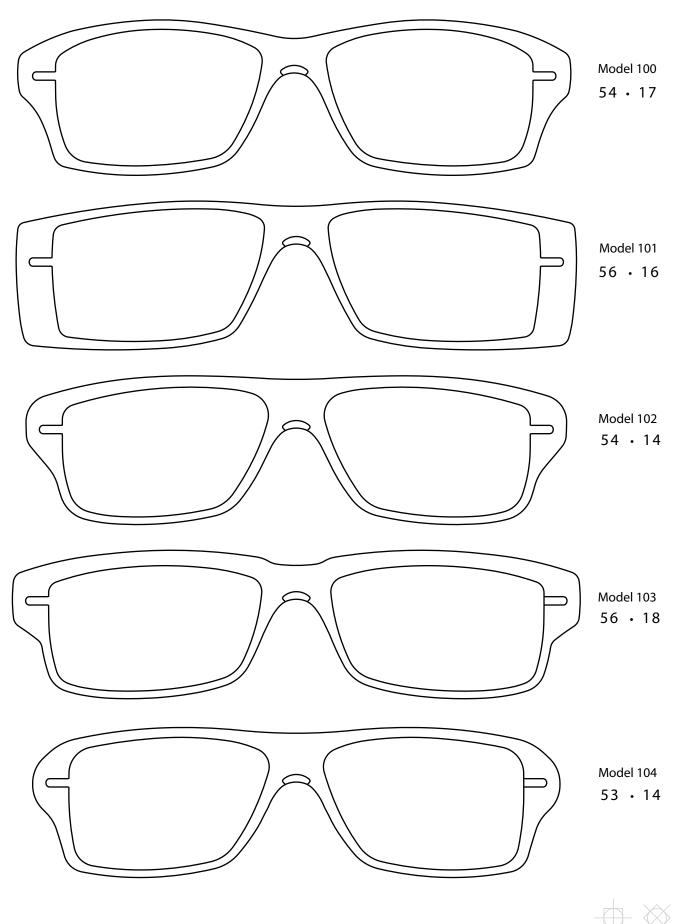



¹⁰⁰V1814

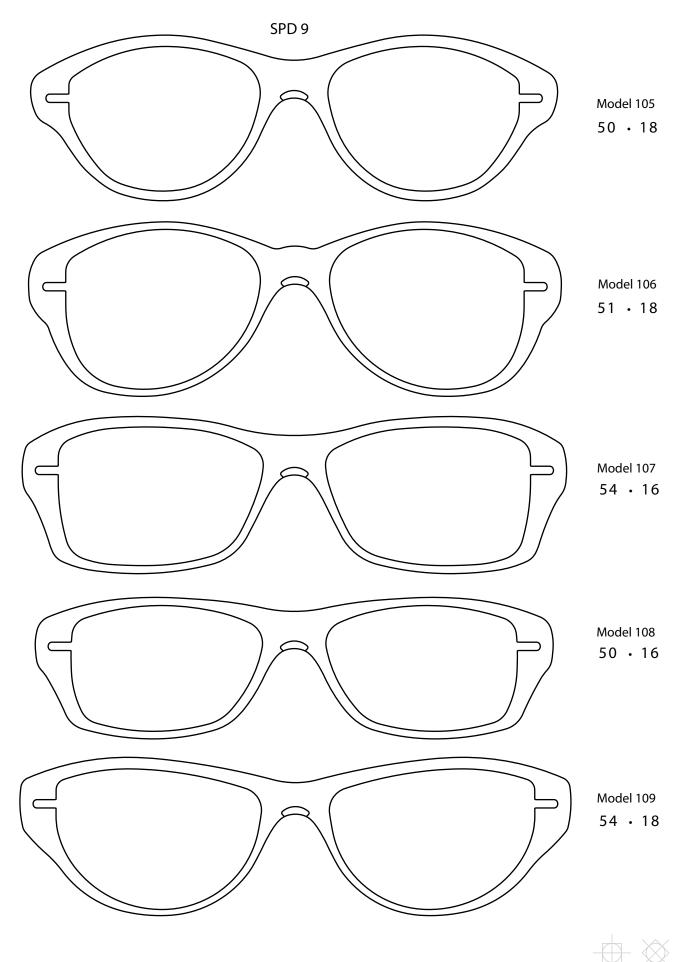

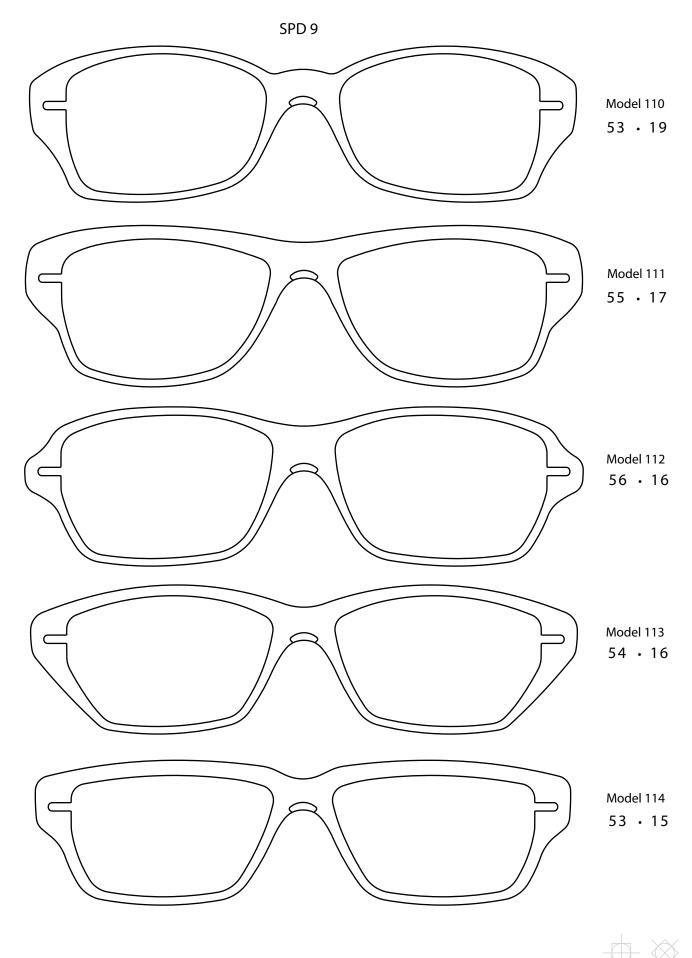

SPD 9

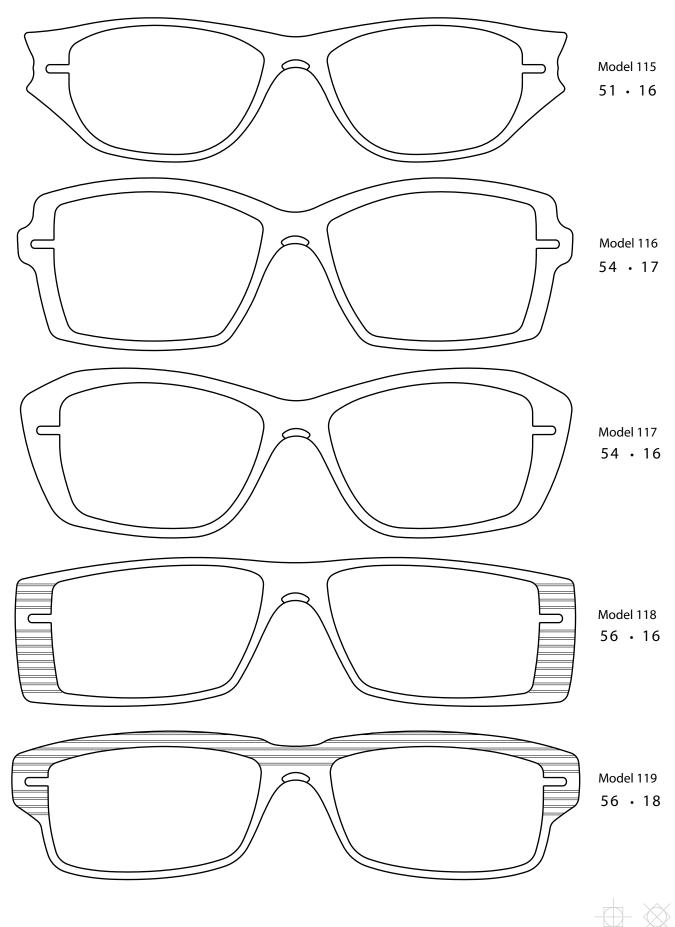


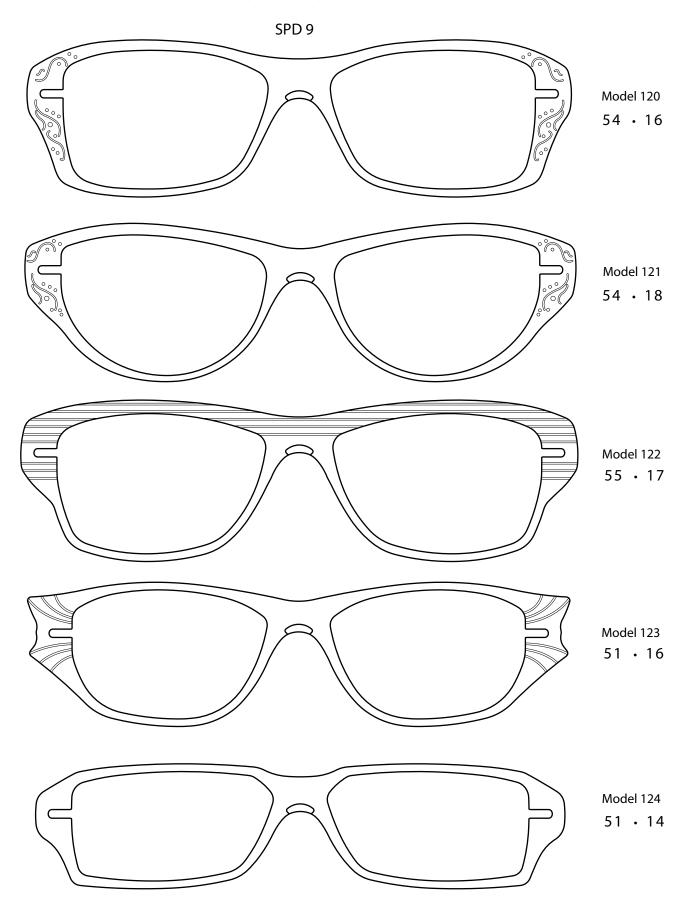




www.polariseyewear.com

SPD 9

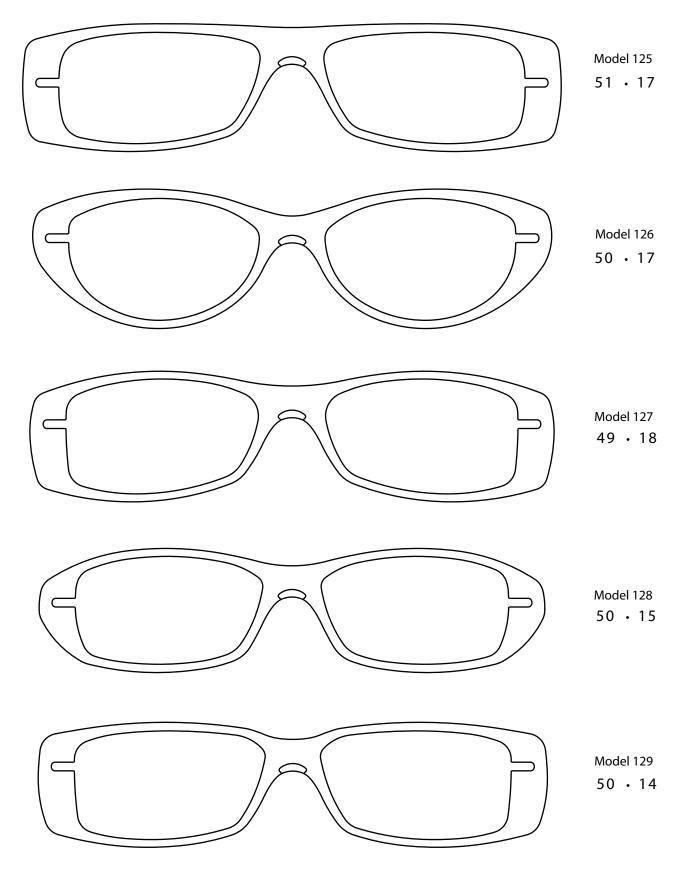



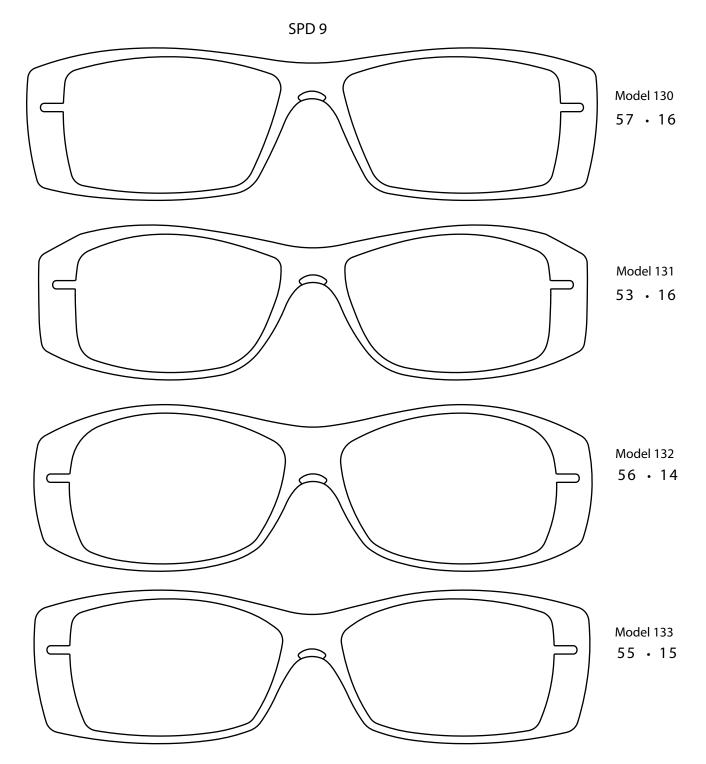



SPD 6






A smooth design in a variety of colours to match your style.

SPD 6 SIDE COLOURS




NOTE: Due to the printing process, colours reproduced in this catalogue may vary from the actual product. Products are not shown to actual size.


FRONT PIECE COLOURS

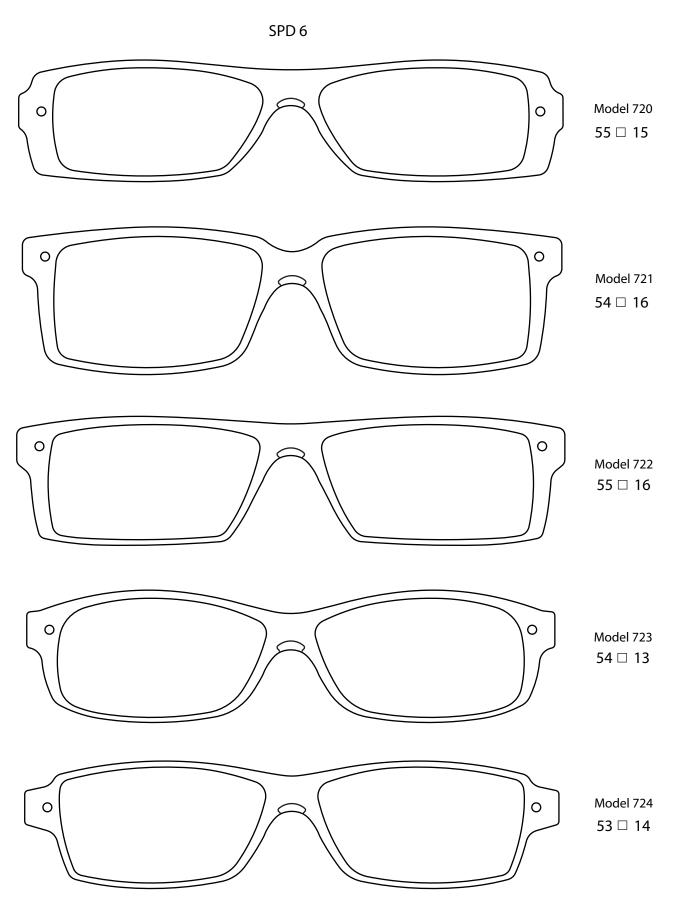



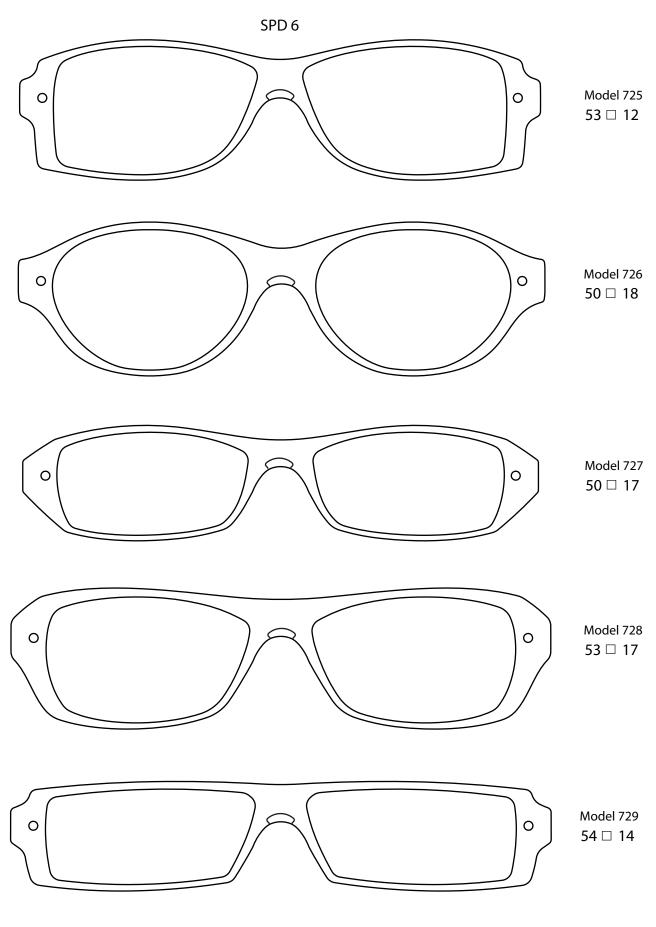

www.polariseyewear.com

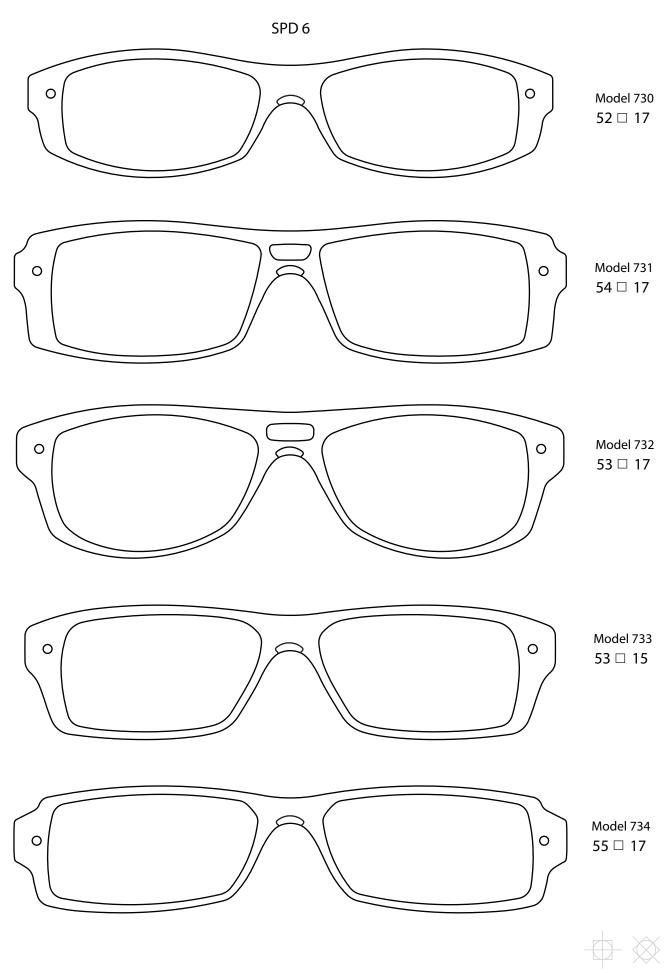

720V5013


SPD FRONT COLOURS



NOTE: Due to the printing process, colours reproduced in this catalogue may vary from the actual product. Products are not shown to actual size.

PATTERN DESIGNS


Patterns designs available on **any shape** in matt and solid colours.

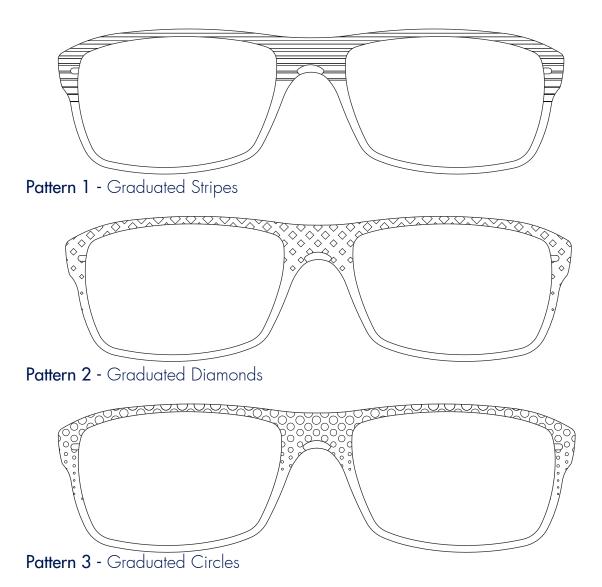


SIZE ALTERATIONS

Staffan Preutz Design can be tailored to fit perfectly for the best cosmetic and optical effect. The size modifications available are listed below:

- ***1** = -2mm each lens
- *2 = -4mm each lens
- ***3** = +2mm each lens
- *4 = +4mm each lens
- ***5** = -2mm in width each lens
- ***6** = -4mm in width each lens
- ***7** = +2mm in width each lens
- ***8** = +4mm in width each lens
- *9 = +2mm on bottom each lens

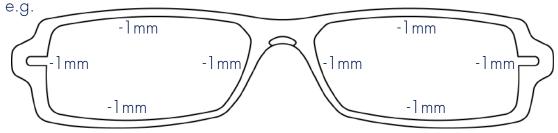

To reduce each lens size by 2mm in width and height as above, order size change *1 (-2mm).

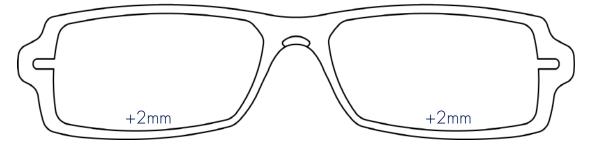

To add 2mm in depth to the shape, useful for accommodating progressive lenses order size change *9 (+2mm on bottom).

NOTE: Distance between lenses (DBL) remains constant when changing size.

MINOR ADJUSTMENTS

- The bridge and lugs can be manipulated using light pressure. No heating is required.

- Length to bend behind the ear can be specified when ordering or easily adjusted in store.

- An alternative bridge is available for those that require the frame to sit slightly higher and further away from the face. For more information please contact our customer service team.

GLAZING INFORMATION

- All lens materials can be used to glaze Staffan Preutz Design with the exception of glass.

- Make a groove around the edge of each lens 0.8mm deep and the same width as a traditional supra frame.

- Assemble each lens by inserting it into the frame from the back, starting with the sharpest corner.

- No heat is required to glaze or adjust Staffan Preutz Design frames.

- For further glazing advice please contact our customer service team.

SAMPLE CODES

Engraved on the lens of a sample is a code referring to the products specifications.

- E.g. SPD P100*1/26/30

- P100 indicates the shape number. The *1 indicates a -2 size. (If standard '0' size, no size code will be displayed after the shape number)

- 26 represents the colour and pattern of the front sheet

- 30 signifies the colour and design of the sides